



# Linux Installation

This guide is specific to the installation of Partek Genomics Suite software on a Linux operating system.


1. Download the latest Linux version of [Partek Genomics Suite](#).
2. Click on Download under the Linux session and save/open the file. Extract on you preferred folder and remember its location.

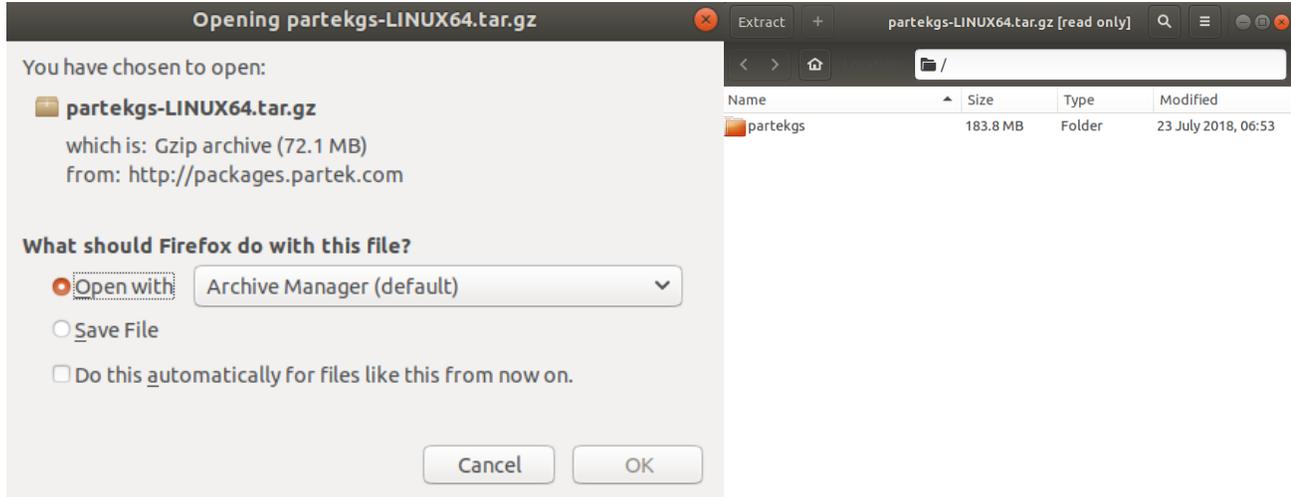

4. Install FlexNet following this [guide](#).
5. Open a terminal and navigate to the location where you extracted the "partekgs" folder.
6. Go into the "bin" directory and run the "partek" file (Figure 2).

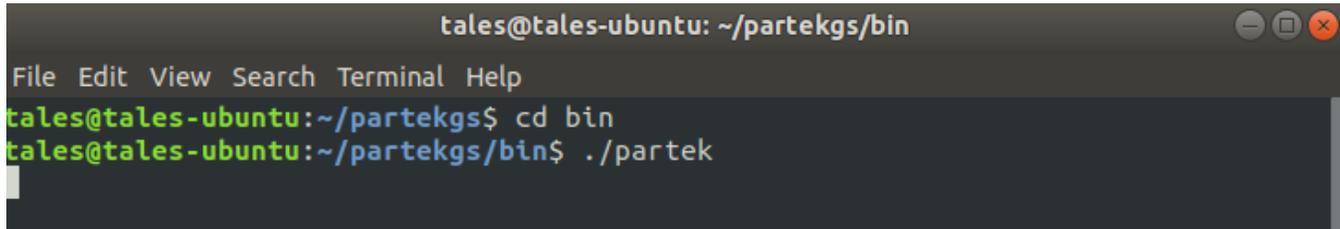

7. The Partek License Manager will appear.

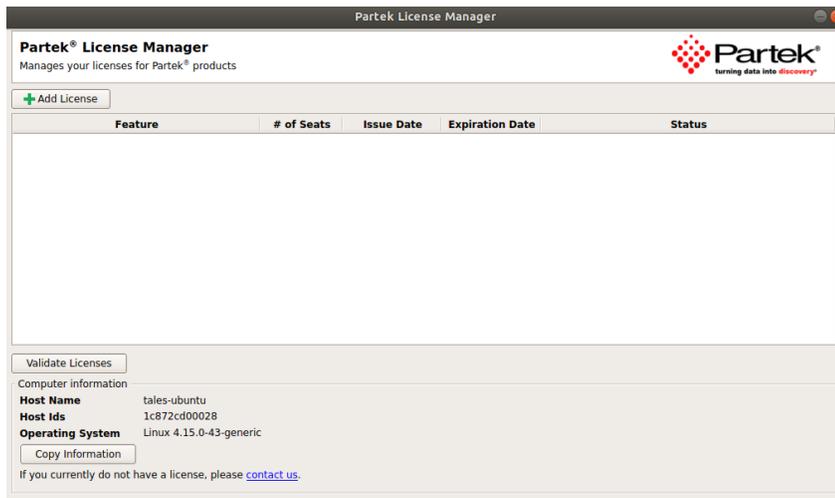

8. If you do not have a license, please contact [licensing@partek.com](mailto:licensing@partek.com).
  - Be sure to click on the "Copy Information" button under the "Computer information" section of this window to paste the information into the email to [licensing@partek.com](mailto:licensing@partek.com).
9. If you have a license, click on the "+ Add License" button.

10. The Partek License Manager - Add License window will appear.

11. Select the "License server" radio button and fill the "Server Name" with the "Host Name" that you copied from the "Computer information" section.

12. Select "Add".

11. After adding the license, the Partek License Manager window will reappear. Select "Validate Licenses" and the license information will show up on your window.




The screenshot shows the Partek License Manager window. At the top, it says "Partek License Manager" and "Manages your licenses for Partek products". There is a logo for Partek with the tagline "turning data into discovery". Below that is a "+ Add License" button. The main area contains a table with the following data:

Feature	# of Seats	Issue Date	Expiration Date	Status
Partek Genomics Suite	1	2019-Feb-06	2019-Dec-31	Valid
Pathway	1	2019-Feb-06	2019-Dec-31	Valid

12. Close the Partek License Manager and run "./partek" again from the "partekgs/bin" directory.

13. The Partek License Agreement window will appear, select "Agree".

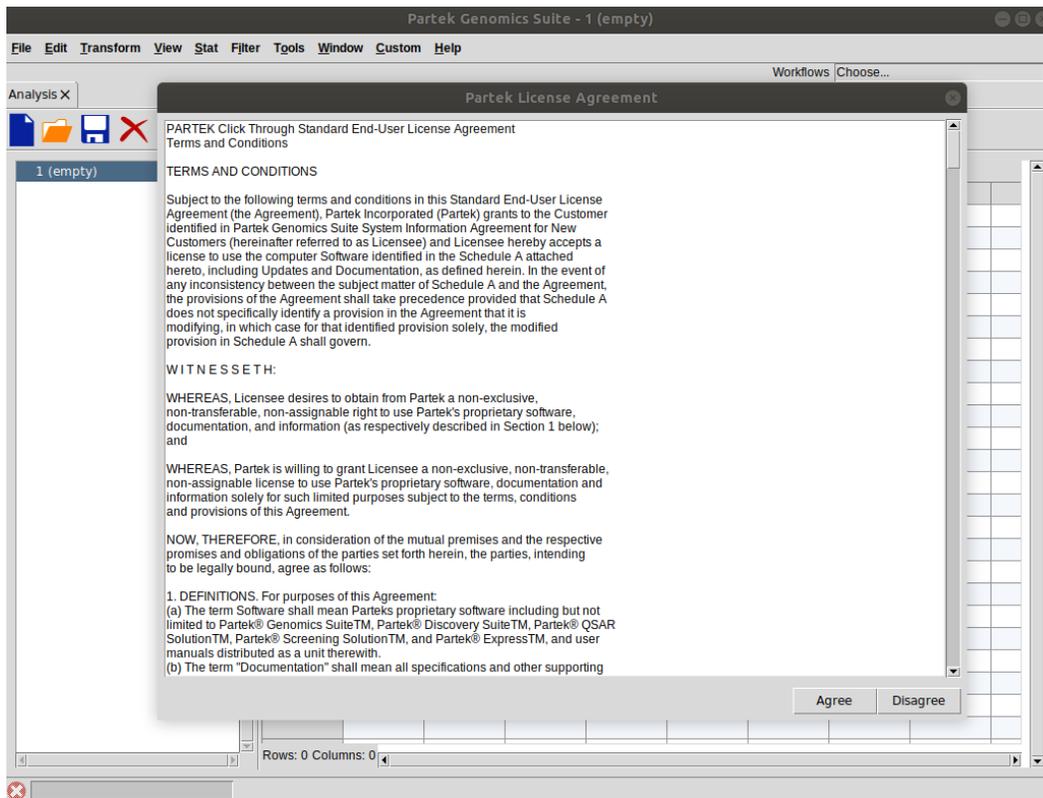

You are now able to use Partek Genomics Suite.

## Additional Assistance

If you need additional assistance, please visit [our support page](#) to submit a help ticket or find phone numbers for regional support.