

Filter reads

Next generation sequencing (NGS) data is notably huge in file size. Dealing with NGS data is not only time consuming but also puts constraints on hard disk space. This is especially true if analysis parameters need to be optimized. The *Filter reads* task is a very useful tool to get a subset of the raw data upon which optimization can be performed. The optimized parameters can then be saved and applied to the whole dataset.

Filter reads is only available for unaligned reads of FASTQ format. Select the *Unaligned Reads* data node then select **Filter reads** from the *Pre-alignment tools* section on the menu.

There are two options to filter reads: **Subsample reads** and **Filter by read length**.

To *Subsample reads*, specify how many reads you want to keep for every nth reads. For example: if the user specifies to "Keep one read for every 10 reads" (Figure 1), this means that for every 10 reads, the program will keep only 1 read. This is equivalent to keeping 10% of the data.

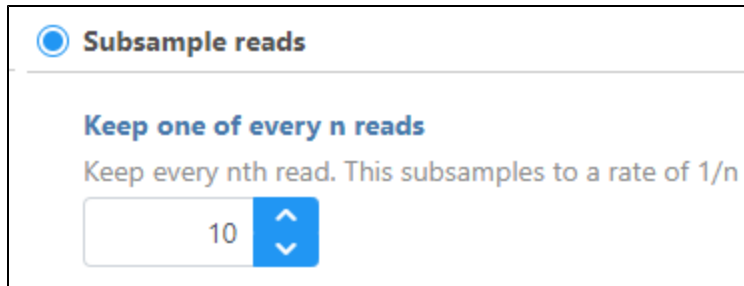
The screenshot shows a web interface for the 'Subsample reads' filter. At the top, there is a radio button icon followed by the text 'Subsample reads'. Below this, the instruction 'Keep one of every n reads' is displayed in blue. Underneath, a greyed-out text says 'Keep every nth read. This subsamples to a rate of 1/n'. A text input field contains the number '10', and to its right is a blue button with up and down arrow icons.

Figure 1. By default, subsample raw data by keeping one read for every 10 reads

To *Filter by read length*, set the read length limits by choosing the minimum and maximum read length(s) to keep.

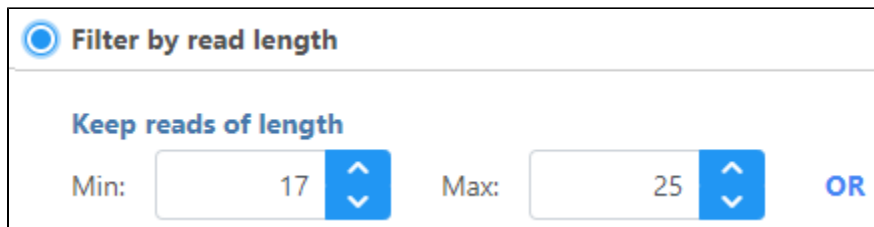
The screenshot shows a web interface for the 'Filter by read length' filter. At the top, there is a radio button icon followed by the text 'Filter by read length'. Below this, the instruction 'Keep reads of length' is displayed in blue. Underneath, there are two input fields. The first is labeled 'Min:' and contains the number '17', followed by a blue button with up and down arrow icons. The second is labeled 'Max:' and contains the number '25', followed by a blue button with up and down arrow icons. To the right of these fields is the word 'OR' in blue.

Figure 2. Filter by read length(s) by setting the parameters for minimum and maximum read length

Additional Assistance

If you need additional assistance, please visit [our support page](#) to submit a help ticket or find phone numbers for regional support.



Your Rating:      Results:      44 rates