

Importing Data

Partek Flow can import a wide variety of data types including raw counts, matrices, microarray, variant call files as well as unaligned and aligned NGS data.

- [Navigating the file browser to transfer files to the server](#)
- [Associate fastq files for multi-omic data](#)
- [SFTP File Transfer Instructions](#)
- [Sample Table from a Text File](#)
- [Import single cell data](#)
- [Importing 10x Genomics Matrix Files](#)
- [Importing and Demultiplexing Illumina BCL Files](#)
- [Partek Flow Uploader for Ion Torrent](#)
- [Importing 10x Genomics .bcl Files](#)
- [Import a GEO / ENA project](#)

The following file types are valid and will be recognized by the Partek Flow file browser.

• bam • bcf • bcl • bgx • bpm • cbcl • CEL • csv • fa • fasta • fastq	• fcs • fna • fq • gz • h5ad • h5 matrix • idat • loom • mtz • probe_tab • qual • raw • rds	• sam • sff • sra • tar • tsv • txt • vcf • zip
---	---	---

In cases where paired end fastq data is present, files will also be automatically recognized and their paired relationship will be maintained throughout the analysis.

Matching on paired end files is based on file names: every character in both file names must match, except for the section that determines whether a file is the first or the second file. For instance, if the first file contains "_R1", "_1", "_F3", "_F5" in the file name, the second file must contain something in the lines with the following: "_R2", "_2", "_F5", "_F5-P2", "_F5-BC", "_R3", "_R5" etc. The identifying section must be separated from the rest of the filename with underscores or dots. If two conflicting identifiers are present then the file is treated as single end. For example, s_1_1 matches s_1_2, as described above. However, s_2_1 does not mate with s_1_2 and the files will be treated as two single-end files.

Apart from paired-end data, files with conventional filename suffixes that indicate that they belong to the same sample are consolidated. These suffixes include:

- Adapter sequences
 - "_bbbbbb" followed by "_" or at the end of the file name, where each "b" is "A", "C", "G", or "T"
- Lane numbers
 - "_L####" followed by "_" or at the end of the file name, where each "#" is a digit 0 to 9
- Dates
 - in the form "####-##-##" preceded or followed by a period or underscore
- Set number
 - of the form "_####" from the end

Navigating the file browser to transfer files to the server


The file browser is used to transfer files to the server so that these files can be added to a project for analysis. If you are importing a Bioproject from GEO /ENA or using URLs for data import, there is no need to transfer the files to the server.

To access the file browser and upload data to the server, use any of these options:

- access **Transfer files**  **Transfer files** on the Partek® Flow® homepage

- within a project, after selecting the file type to transfer, using the **transfer files link** [Transfer files to the server](#) available within all file import options
- from the settings, go to **Access management > Transfer files**

Using the file browser to transfer files to the server:

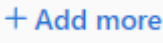- Click **Transfer files**  to access the file browser
- Drag and drop or click **My Device**  to add files from your machine
- Click **Browse**  to modify the *Upload directory* or create a new folder. The **Upload directory**  should be specified, known, and distinguishable for project file management. You will return to this directory and access the files to import them into a project
- To continue to add more files use **+ Add more**  in the top right corner. To cancel the process select **Cancel** in the top left corner
- Click **Upload**  to complete the file upload
- Do not exit the browser tab or let the computer go to sleep or shut down until the transfer has completed

File size displayed in the table is binary format, not decimal format (e.g. GB displayed in the table is gigabyte not gigabyte. 1 gigabyte is 1,073,741,824 bytes. 1 gigabyte is 1,000,000,000 bytes. 1 gigabyte is 1.074 gigabytes).

Associate fastq files for multi-omic data

There are projects with more than one file type, such as single cell multi-omic assays that generate protein and RNA together. In these cases, files need to be associated with each other if starting in fastq format. If we start with processed data, there is no need to associate these files.


- Define the type of data the file represents when importing files into the project.

☒ **Select data type: (only required for multi assay project)**

Select your data type. If none of the preset options matches your data, choose custom and add a custom type.

ATAC-Seq ▼

If this step is skipped, the data type can be changed after import by right clicking the data node.

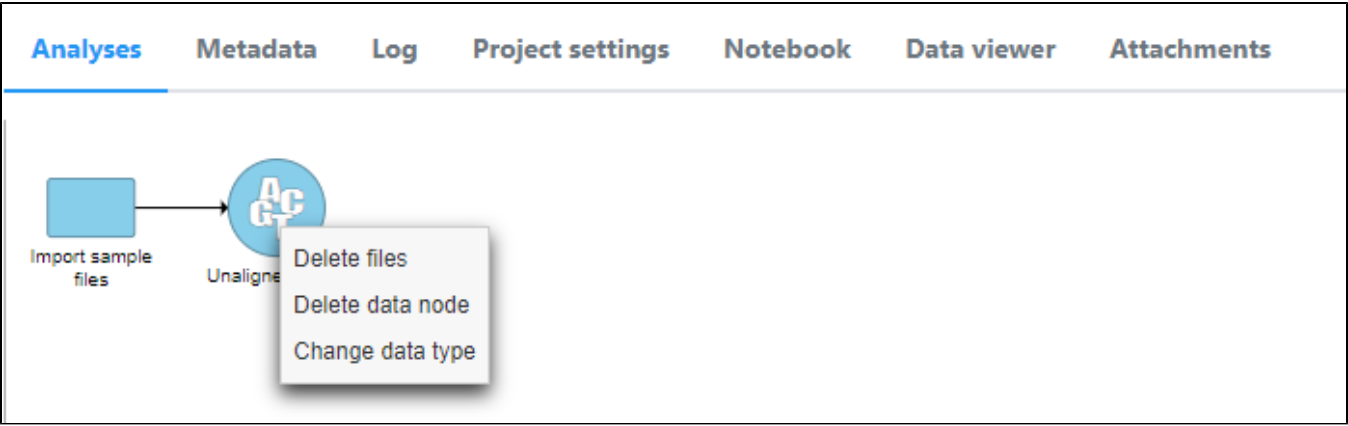

- After importing both types of data, associate fastqs with the already imported data. (e.g. associate RNA fastqs with ATAC fastqs already imported into Partek Flow.

Analyses	Data	Log	Project settings	Notebook	Data viewer	Attachments
Sample attributes		Sample name		Attributes		Files
Manage				Type	Subject ID	fastq
Assign values		1	ERR1618872	Tumor	ESP2	2 files: R136360_ERR1618872_1, R136361_ERR1618872_2
		2	ERR1619078	Normal	ESP2	2 files: R136363_ERR1619078_1, R136364_ERR1619078_2
		Associate files with this sample				

Additional Assistance

If you need additional assistance, please visit [our support page](#) to submit a help ticket or find phone numbers for regional support.

★

★

★

★

★

★

★

★

★

★

Your Rating: ★★★★★

Results: ★★★★★ 35 rates

Partek